
The EasyLED series integrates advanced LED illumination and control functions into a single compact device. It is the ideal standard illumination system for routine examinations, research, and education applications.



### Spot illumination

The EasyLED Spotlight and the EasyLED Double Spotlight Plus are an ingenious design that allows mounting to any microscope stand.

The intensity control is ergonomically located near the microscope focus knob. The robust metal lighthead emits up to 130 lm of cool light and rated lifetime is 50,000 hours minimum.

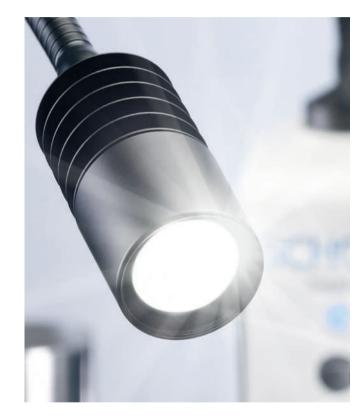


#### Transmitted illumination

The EasyLED Transmitted Light Stage features a SCHOTT Opalika® surface that offers homogeneous and shadow free illumination.

Its robust metal housing and a functional design, delivers up to 12,000 cd/m<sup>2</sup> and a long life time of at least 50,000 hours



SCHOTT Authorised Distributor
Lighting and Imaging Products
Silicon Connection Pte Ltd
Blk 7 Kallang Place
#05-01/05/07/11/12
Singapore 339153
Phone +65 81861235
steven@si-cnx.com
www.si-cnx.com/products/lighting-imaging/

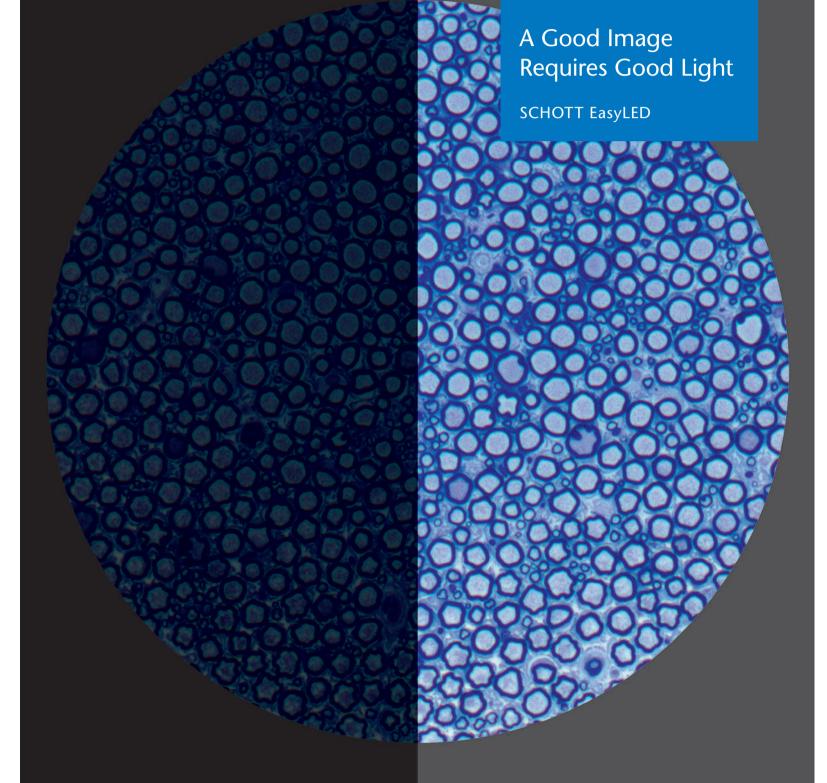

# A good image requires good light

#### SCHOTT EasyLED illumination system

Optimal illumination of objects under examination is one of the most decisive factors for good results in stereo-microscopy. Ideal light is necessary for detection of hidden details, differentiation of specific areas and emphasis of special characteristics.





#### Ring illumination

The EasyLED Ringlight enables homogeneous and shadowless illumination with a maximum brightness of 90 kLux and a long life of at least 50,000 hours. It features a robust metal housing and a highly functional design.

The powerful 140 kLux EasyLED Ringlight Plus has additional controllable segments that enable new contrast methods. The integrated control element allows easy and intuitive switching between the different illumination modes, as well as a rotation function in both directions.



We offer suitable accessories for the EasyLED Ringlight Plus for the reduction of disturbing surface reflections.

The ringlight diffuser enables homogenization and diffusion of the light, which significantly reduces surface reflections. The Polarizer Set, which consists of a rotating polarizer in front of the LEDs and an analyzer in the observation ray path, makes it possible not only to eliminate reflections, but also to emphasise double-refracting materials.

Caption title: Cross section through the sciatic nerve, illuminated with EasyLED Transmitted Light Stage.



FSC\* C006655

30074 ENGLISCH 09/2019 Ali/odd Printed in Germany







Oblique lighting


LED Ringlight and polarization

Oblique lighting





LED Ringlight

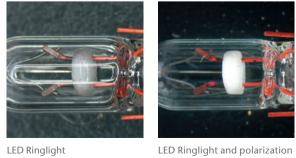

## EasyLED Ringlight

- P/N 600.200
- 45 high-power LEDs • Working distance 50-135 mm
- Integrated control
- Lifetime > 50,000 hours

High intensity (90 kLux)

## EasyLED Ringlight Plus

P/N 600.300

- Highest intensity (140 kLux)
- 48 high-power LEDs
- Working distance 50-130 mm
- Integrated control Segmentable
- Lifetime > 50,000 hours

#### Accessories

- For EasyLED Ringlights
- To reduce surface reflections


Diffuser

P/N 400.570

Polarizer Set

P/N 400.550



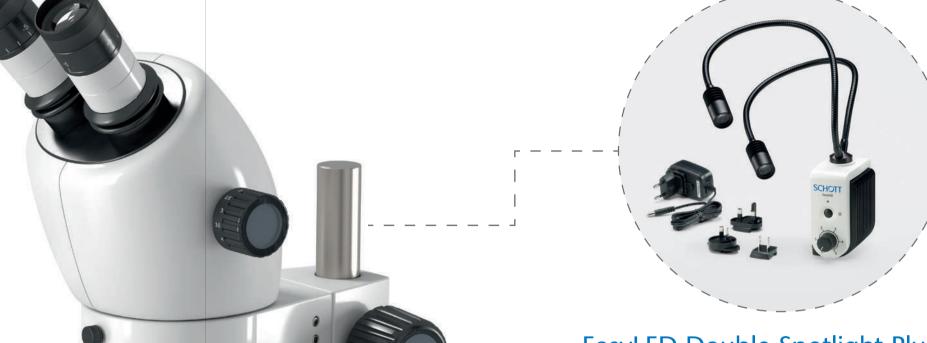

EasyLED Double Spotlight Plus





### EasyLED Double Spotlight Plus P/N 600.121

• Intensive 130 lm per lighthead

EasyLED Transmitted Light Stage


• High luminance of 12,000 cd/m<sup>2</sup>

P/N 600.400

Integrated control

 Scratch-resistant surface • Lifetime > 50,000 hours

- Lightheads can be controlled individually
- Lifetime > 50,000 hours

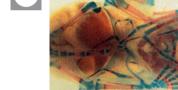

LED Ringlight






LED Ringlight and diffuser


Segmentation LED Ringlight





Segmentation LED Ringlight

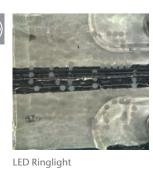

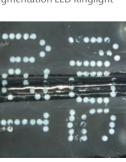












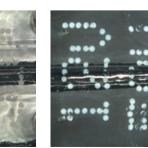

LED Ringlight and polarization



